
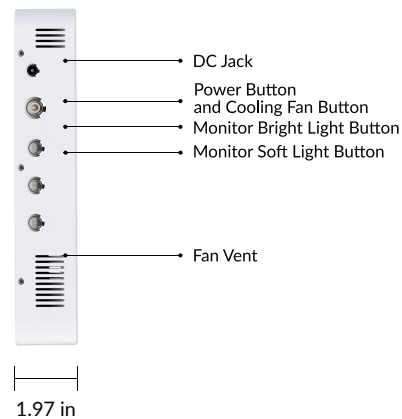
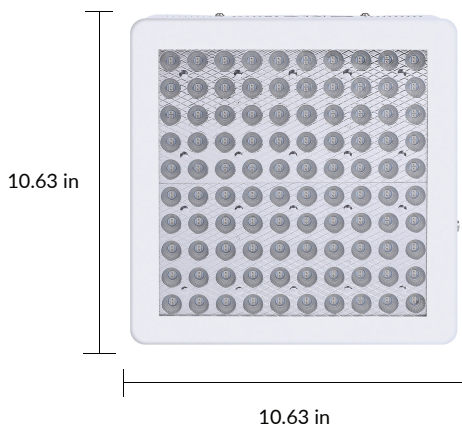
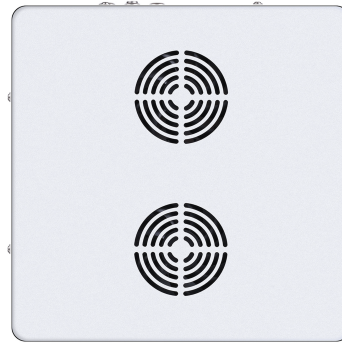


PP100

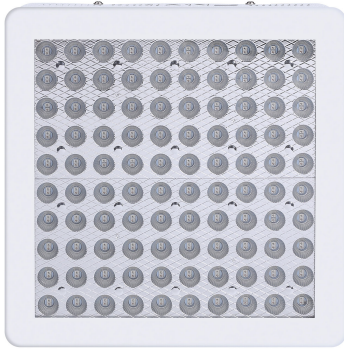

LED Illuminator Specification:

LED Power Class:	100w
Dimensions:	10.63"x10.63"x1.97"
Power Input:	DC 24V
Power Consumption:	85-100w
LEDs:	200x0.5w
Focusing Lens Bean Angle:	60-90 Degree
Cooling Fans:	2 Fans
Lifespan:	100,000 Hrs
EMF Emission:	0.0 μ T
Light Output:	70-100mw/cm ²

User Guide

1. Plug DC power adapter into the DC Jack on the device.
2. Push power button (green) to turn on unit.
3. Use 850nm wavelength, push blue button.
4. Use 660nm wavelength, push red button.
5. Use 40HZ frequency button to active frequency model.

Feature

- 24V DC Power Source Input
- Easy to carry for safety concern
- 100W Compact Size Version
- 200pc 2835 LED Chip Special Design
- Power Output is distributed balance

■ Connection:

The light therapy device has a wired power input socket.
The center of the socket is positive pole.
The socket is connected to 24V DC / AC power supplier.

FCC NOTICE

This equipment has been tested and found to comply with the limits for a class of a digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates radio frequency, and if not installed and used in accordance with the instruction manual, may cause harmful interference.

Warranty

PROFESSIONAL COLLECTION LED Technologies, warrants that this product will be free from defects in materials and workmanship for 1 year from the date of purchase by the original owner. This warranty covers the original purchaser only and is not assignable or transferable to another party. Register your warranty online at <https://cmvisionLightTherapy.com> or by calling CustomerService at 713-779-8788.



The lightning flash with arrowhead symbol, with an equilateral triangle, is intended to alert the user to the presence of "uninsulated" dangerous voltage" within the product's enclosure that may be of sufficient magnitude to constitute a risk of electric shock to persons.



This exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the appliance.



CE Mark
This apparatus is manufactured to comply with the radio interference.

Benefits

- Anti-Aging
- Fat Loss
- Muscle Gain
- Performance/Recovery
- Brain Optimization

■ Notes:

To ensure longer lifetime and more efficient, this led illuminator must be used with stable 24V DC / AC power source. The voltage range should be within 24V DC/AC +5%, if voltage is too high the illuminator lights will be damaged, if voltage is too low, the illuminator lights will not work properly.

Customer Service

You can reach LED Technologies Customer Service by calling **713-779-8788**
Monday-Friday 8:00am-5:00pm,
Mountain Standard Time.
Customer Service can also be reached by email at support@cmvisionlighttherapy.com
Our Website: <https://cmvisionLightTherapy.com>

Disclaimer

This product is not intended to diagnose, cure or prevent disease. LED Technologies makes no claims, representations or warranties regarding the ability of this product to cure any physical, skin or mental condition. A qualified health professional should be consulted with regard to any condition requiring medical attention.



ROHS Announcement (Option)
All lead-free products offered by the company comply with the requirements of the European law on the Restriction of Hazardous Substances (RoHS) directive, which means our manufacture processes and products are strictly "lead-free" and without the hazardous substances cited in the directive.



The crossed-out wheeled bin mark symbolizes that within the European Union the product must be collected separately at the product end-of-life. This applies to your product and any peripherals marked with this symbol. Do not dispose of these products as unsorted municipal waste.

The company does not warrant that this manual will be uninterrupted or error-free.
We reserve the right to revise or remove any content in this manual at any time.